



Shoulder Wrap



NIPPON SIGMAX CO.,LTD.
33rd Floor Shinjuku Oak Tower
6-8-1 Nishi-shinjuku, Shinjuku-ku,
Tokyo 163-6033 JAPAN
Phone: +81 3 5326 3229

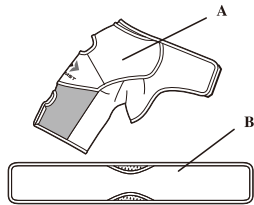


1

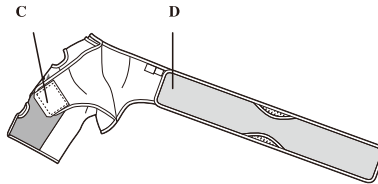


2

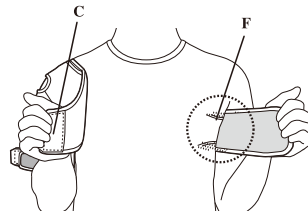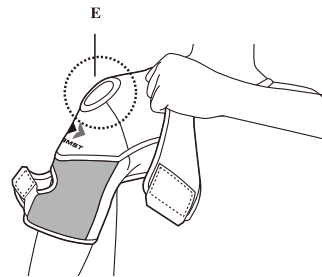

4

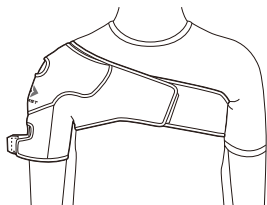

5

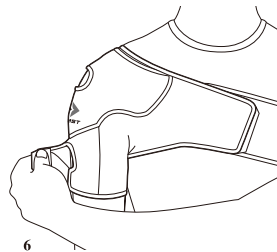

6

ZAMST SHOULDER WRAP INSTRUCTIONS FOR USE

Please read through the following instructions completely before using the product.

Intended use: Shoulder compression and protection

*This product is not for medical treatment.

*This product is for shoulder support; however, it will not prevent all shoulder injuries.

Do not use this product in the following cases:

- If you have abnormal conditions such as cuts, fractures, numbness, swelling or rash in the area where product is used.

Safety precautions:

- If you experience pain, swelling, skin irritation, discomfort or any other abnormal reaction while using this product, discontinue use and consult with a physician.

- If your skin is sensitive or you have allergies, you may develop a rash or allergic reactions while using the product. If such symptoms occur, discontinue use and consult with a physician.

- Do not use the product while sleeping. If the product shifts to an abnormal position or is used for extended periods while sleeping, it may interfere with blood circulation.

- If the body part is compressed for extended periods, blood circulation may be adversely affected or you may experience discomfort. To avoid such conditions, it is recommended that you remove and re-apply the product at regular intervals.

- If the product is broken or worn or if the hook-and-loop fastener does not work, do not use the product.

- Make sure that you use the correct size. (Refer to the size chart on the package and size tag inside the product.)

- Do not pull the strap too tight or create excessive compression with the support.

- If the product shifts out of place or becomes too loose, re-apply the product starting with the first step described in the "How to apply" instructions.

- Do not modify or repair the product. The product may lose its effectiveness and/or cause injuries.

To prevent product damage:

- The product may become discolored by perspiration or other liquids. The product's dye may bleed onto clothing worn under or near the product.

- Make sure that the hook-and-loop fasteners do not touch fabric. They may damage fabric.

How to apply:

*Illustrations are for the right shoulder. The product is not designed specifically for use on the right or left, so application is the same for the left shoulder.

1. The product consists of a shoulder wrap (A) and a chest strap (B). Release all hook-and-loop fasteners. (Illustration 1)

2. The shoulder wrap has hook fasteners on the chest side (C) and the back side (D). Attach the chest strap to the hook on the back side of the shoulder wrap. (Illustration 2)

3. Slide the shoulder wrap up onto the arm and shoulder until the point of the shoulder is positioned under the hole (E). (Illustration 3)

4. Pull the chest strap from behind and pass it under the opposite shoulder (F). (Illustration 4)

5. Pull the chest strap across the chest with the appropriate compression. Fasten the hook at the chest. (Illustration 5)

6. Pull the strap on the upper arm with appropriate compression. Fasten hook. (Illustration 6)

How to wash:

- Hand wash with mild detergent and air dry.

- Do not place in dryer. Do not dry clean or iron. Do not use products containing chlorine, including bleach. These may damage the product.

- The product's dye may bleed. Do not wash with whites or bright colors.

Composition: Nylon, Polystyrene, Polyester, Polyurethane, Cotton